Homedics

Foot flow


✓ Register your product today at
www.homedics.co.uk/product-registration

3 YEAR GUARANTEE


PSL-2000H-EB

PRODUCT FEATURES

1. Remote Control
2. Vibrating Foot Plates
3. Digital Display
4. Mains adapter with interchangeable pins
5. Jack socket

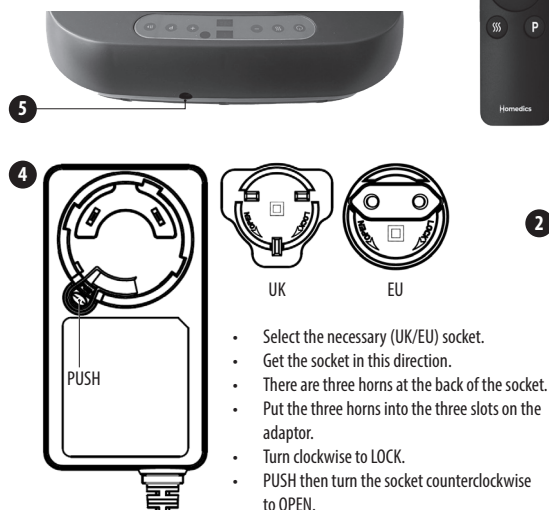- Select the necessary (UK/EU) socket.
- Get the socket in this direction.
- There are three horns at the back of the socket.
- Put the three horns into the three slots on the adaptor.
- Turn clockwise to LOCK.
- PUSH then turn the socket counterclockwise to OPEN.

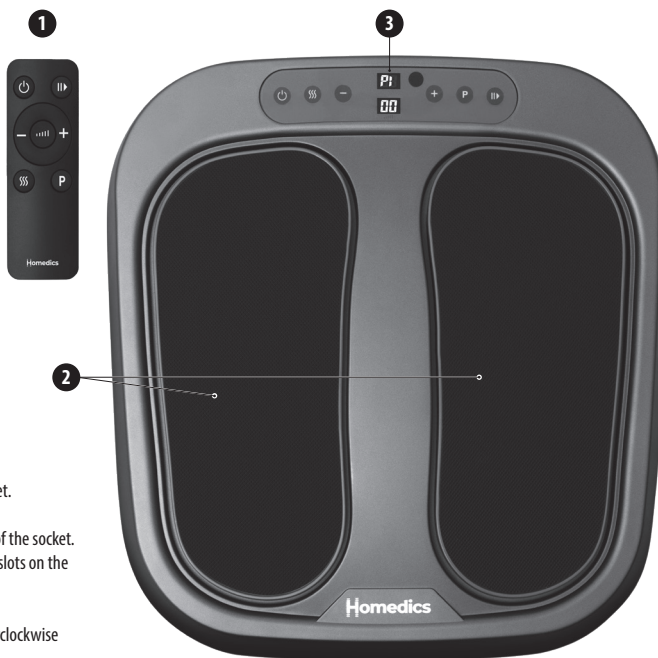CONTROL BUTTONS AND DIGITAL DISPLAY

Control Buttons

- Power** – Press to switch the device on and off
- Program** – Press once to select **Manual Mode P1**, and again to select **Program P2** or **P3**
- Pause / Play** – Press to pause or restart the current program
- Minus** – Press to decrease speed (works in Manual Mode only)
- Plus** – Press to increase speed (works in Manual Mode only)
- Heat** – Press twice for **low heat** (level 1), and again for level 2

Understanding the Digital Display

P1 Manual Mode P1

06 Vibration Speed 6

NOTE: In Manual Mode you can adjust the vibration speed with the + and - buttons

14 14 mins remaining

09 Vibration Speed 9

NOTE: The entire display will flash when the device is paused

HE Heat
00 OFF

HE Heat
01 level 1

P2 Program P2

18 Vibration Speed 18

NOTE: When program P2 or P3 is selected, the vibration speed will vary automatically throughout the program

HE Heat
02 level 2

BEFORE FIRST USE

1. Remove all packaging
2. Remove the plastic protective strip from the remote control to activate the battery
3. Slide the appropriate pins fully into the mains adapter plug, either 3-pin GB, or 2-pin EU.
4. Push the jack plug end of the cable into the jack socket at the rear of the device.
5. Place Foot Flow on the floor in front of your chair, ensuring the surface is flat and level
6. Plug the mains adapter into a nearby mains socket.

INSTRUCTIONS FOR USE

1. Sit in a chair and place your feet on the foot plates. Never stand on the footplates, or use the device on one side only.
2. Press the power button to switch Foot Flow on.
NOTE: The device will start up in manual mode (P1) and on vibration speed 6
3. Adjust the vibration speed as required with the + and - buttons
4. Add heat, if desired, by pressing the heat button twice for low heat, level 1, press again for level 2. Allow a minute or two for the foot plates to warm up fully.
5. Sit back and relax. The device will automatically switch off after 15 mins. In off mode, this product will consume power at a rate of 0.2W.

6. The device can be paused (with the pause/play button) or switched off (with the power button) at any time.

NOTE: All functions can be controlled either with the remote control or the buttons on the device. Recommended use time is 15 mins.

Using the programs

To select a program, first switch the device on, then press the Program button until either P2 or P3 appears on the display. Both programs are 15 minutes long and will automatically vary the intensity settings throughout the program for an invigorating treatment.

Changing the remote control battery

1. Remove the screw on the back of the remote control
2. Open the cover and remove the battery
3. Insert a new CR2032 battery with the + marking facing upwards
4. Close the cover and replace the screw.

Cleaning your device

- Make sure the device is unplugged from any power supply and allow to cool before cleaning. Clean only with a soft, slightly damp sponge.
- Never allow water or any other liquid to come into contact with the appliance.
- Do not immerse in any liquid to clean.
- Never use abrasive cleaners, brushes, glass/furniture polish, paint thinner, etc to clean.

Storage

Place the appliance in its box or in a safe, dry, cool place. To avoid breakage, DO NOT wrap the power cord around the appliance. DO NOT hang the unit by the cord.

READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- NEVER leave an appliance unattended when plugged in. Unplug from the mains outlet when not in use and before putting on or taking off parts or attachments.
- This appliance is suitable for ages 8 and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins, metallic fasteners or objects into the appliance or any opening.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by HoMedics.
- NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the HoMedics Service Centre for examination and repair.
- Excessive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- If the product feels excessively hot, switch off at mains and contact the HoMedics Service Centre.
- Keep cord away from heated surfaces.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the cord or use cord as a handle.
- Do not run cord under carpeting. Do not cover cord with throw rugs, runners or similar coverings. Do not route cord under furniture or appliances. Arrange cord away from traffic areas where it can be tripped over.
- DO NOT use outdoors.
- This product requires a 220-240 V AC power supply.
- The appliance is only to be used with the power supply unit provided with the appliance.
- DO NOT attempt to repair the appliance. There are no user serviceable parts. All servicing of this appliance must be performed at an authorised HoMedics Service Centre.
- DO NOT sit or stand on the footplates. The unit

should only be used when placed on the floor and user is seated.

- The total weight limit for this product is 20kg.
- Never block the air openings of the appliance. Keep air openings free of lint, hair etc.
- If you have any concerns regarding your health, consult a doctor before using this appliance.
- Pregnant women, diabetics and individuals with pacemakers should consult a doctor before using this appliance. Not recommended for use by individuals with sensory deficiencies including diabetic neuropathy.
- DO NOT use on an infant, invalid or on a sleeping or unconscious person.
- This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the lower half of their body.
- NEVER use directly on swollen or inflamed areas or skin eruptions.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your GP.
- DO NOT use as a substitute for medical attention.
- Do not use for longer than the recommended time.
- DO NOT use before going to bed. The massage has a stimulating effect and can delay sleep.
- Only gentle force should be exerted against the footplates in order to eliminate risk of injury.
- Please ensure that all hair, clothing and jewellery are kept clear of the footplates at all times.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- Keep children and pets away from the appliance.
- The device is designed for use with feet, never use it on your back, hips, neck or other parts of the body.
- Failure to follow the above may result in the risk of fire or injury.

3 YEAR GUARANTEE

FKA Brands Ltd guarantees this product from defect in material and workmanship for a period of 3 years from the date of purchase, except as noted below. This FKA Brands Ltd product guarantee does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorised accessory; alteration to the product; or any other conditions whatsoever that are beyond the control of FKA Brands Ltd. This guarantee is effective only if the product is purchased and operated in the UK / EU. A product that requires modification or adaptation to enable it to operate in any country other than the country for which it was designed, manufactured, approved and / or authorised, or repair of products damaged by these modifications is not covered under this guarantee. FKA Brands Ltd shall not be responsible for any type of incidental, consequential or special damages.

To obtain guarantee service on your product, return the product post-paid to your local service centre along with your dated sales receipt (as proof of purchase). Upon receipt, FKA Brands Ltd will repair or replace, as appropriate, your product and return it to you, post-paid. Guarantee is solely through HoMedics Service Centre. Service of this product by anyone other than the HoMedics Service Centre voids the guarantee. This guarantee does not affect your statutory rights.

For your local HoMedics Service Centre, go to

www.homedics.co.uk/servicecentres

WEEE explanation

 This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

Battery directive

 This symbol indicates that batteries must not be disposed of in the domestic waste as they contain substances which can be damaging to the environment and health. Please dispose of batteries in designated collection points.

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FKA Brands Ltd, Somerhill Business Park, Tonbridge, Kent TN11 0GP, UK
FKA Brands Ltd, 29 Earlsfort Terrace, Dublin D02AY28, Ireland